EN	ENGLISH	English	Dental Zirconia
FR	FRANÇAIS	French	Zircone dentaire
ES	ESPAÑOL	Spanish	Óxido de zirconio dental
IT	ITALIANO	Italian	Ossido di zirconio dentale
NL	NEDERLANDS	Dutch	Dentaaal zirkoniumoxide
DE	DEUTSCH	German	Dentales Zirkonoxid
SV	SVENSKA	Swedish	Dental zirkoniumdioxid
NO	NORSK	Norwegian	Dental zirkoniumdioksidi
FI	SUOMI	Finnish	Hammaslääketieteellinen zirkoniumdioksidi
DA	DANSK	Danish	Dentalt zirkonia
PT	PORTUGUÊS	Portuguese	Óxido de zircónio dentário
EL	ΕΛΛΗΝΙΚΑ	Greek	Οδοντιατρικό οξείδιο του ζιρκονίου
TR	TÜRKÇE	Turkish	Dental zirkonyum oksit
PL	POLSKI	Polish	Tlenek cyrkonu do aplikacji dentystycznych
RO	ROMÂNĂ	Romanian	Oxid de zirconiu dentar
HR	HRVATSKI	Croatian	Dentalni cirkonij
HU	MAGYAR	Hungarian	Fogászati cirkonoxid
SL	SLOVENŠČINA	Slovenian	Dentalni cirkonij
CS	ČESKY	Czech	Dentální zirkonoxid
SK	SLOVENSKY	Slovak	Dentálny oxid zirkónu
BG	БЪЛГАРСКИ	Bulgarian	Стоматологичен цирконий
UK	УКРАЇНСЬКА	Ukrainian	Стоматологічний цирконій
ET	EESTI	Estonian	Dentaaalne tsirkoonium
LV	LATVIEŠU	Latvian	Zobārstniecības cirkonijs
LT	LIETUVIŠKAI	Lithuanian	Odontologinis cirkonio oksidas

	EN	UK Responsible Person	PL	Osoba odpowiedzialna w Wielkiej Brytanii
	FR	Personne responsable au Royaume-Uni	RO	Persoana responsabilă din Marea Britanie
	ES	Responsable en el Reino Unido	HR	Osoba odgovorna za Ujedinjeno Kraljevstvo (UK)
	IT	Responsabile per il Regno Unito	HU	Felelős személy az Egyesült Királyságban
	NL	Verantwoordelijke in het VK	SL	Odgovorna oseba v Združenem kraljestvu
	DE	Für UK verantwortliche Person	CS	Odpovědná osoba ve Spojeném království
	SV	Ansvarig person i Storbritannien	SK	Zodpovedná osoba v Spojenom kráľovstve
	NO	Ansvarlig person i Storbritannia	BG	Отговорно лице за Обединеното кралство
	FI	Vastuuhenkilö Yhdistyneessä kuningaskunnassa	UK	Відповідальна особа у Великій Британії
	DA	Ansvarlig person i Det Forenede Kongerige	ET	Suurbritannia vastutav isik
	PT	Responsável no Reino Unido	LV	Apvienotās Karalistes atbildīgā persona
	EL	Υπεύθυνο Πρόσωπο για το Ηνωμένο Βασίλειο	LT	Atsakingas asmuo Jungtinėje Karalystėje
	TR	Birleşik Krallık Sorumlusu		

Kuraray Noritake Dental Inc.
300 Higashiyama, Miyoshi-cho, Miyoshi, Aichi 470-0293, Japan
TEL +81-561-32-8953 FAX +81-561-32-8976

Kuraray Europe GmbH (EU Importer)

Philipp-Reis-Str. 4,
65795 Hattersheim am Main, Germany
Phone: +49 (0)69 305 35835 Fax: +49 (0)69 305 98 35835
URL: <https://www.kuraraynoritake.eu>


EC REP

EMERGO EUROPE

Westervoortsedijk 60, 6827 AT Arnhem, The Netherlands

UTML (Ultra Translucent Multi Layered)

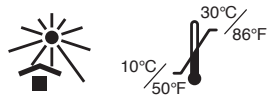I. Introduction

This IFU is for KATANA Zirconia UTML (Ultra Translucent Multi Layered). KATANA Zirconia UTML is a pre-sintered zirconia disc 98.5 mm in diameter, which contains a plastic ring. This is designed for all milling systems using this generic-type disc. (Please refer to your milling system's technical instructions for correct machine operation.) KATANA Zirconia UTML has 2 available thicknesses: (14mm and 18mm) and has 20 shade variations: (A1, A2, A3, A3.5, A4, B1, B2, B3, B4, C1, C2, C3, C4, D2, D3, D4, EA1, EA2, EA3, ENW). KATANA Zirconia UTML consists of 4 gradated shade layers. KATANA Zirconia UTML is recommended for use in fabricating FCZ (Full Contour Zirconia) restorations or the frameworks. The general clinical benefit of this product is to restore tooth function for the following INTENDED USE/ INDICATIONS FOR USE.

II. INTENDED USE/ INDICATIONS FOR USE

KATANA Zirconia is used for the fabrication of the all-ceramic restorations (frameworks, FCZ crowns, FCZ bridges, inlays, onlays and veneers).

[Intended user]

- Dentists (chairside fabrication of restorations; clinical workflow)
- Dental laboratory technicians (fabrication of restorations in the dental laboratory)

III. Sintering Program

Sintering program 1

Temperature	Programming Rate	Holding Time
Room Temp. — 1550°C (2822°F)	10°C/min. (18°F/min.)	—
1550°C (2822°F)	—	2 hrs
1550°C (2822°F) — Room Temp.	-10°C/min. (-18°F/min.)	—

Sintering program 2

Temperature	Programming Rate	Holding Time
Room Temp. — 1400°C (2552°F)	50°C/min. (90°F/min.)	—
1400°C (2552°F) — 1500°C (2732°F)	4°C/min. (7°F/min.)	—
1500°C (2732°F) — 1560°C (2840°F)	10°C/min. (18°F/min.)	—
1560°C (2840°F)	—	16 min
1560°C (2840°F) — Room Temp. (*)	-50°C/min. (-90°F/min.)	—

* The restorations may be removed at 800°C (1472°F) or less depending on the circumstances.

Sintering program 3

Temperature	Programming Rate	Holding Time
Room Temp. — 1450°C (2642°F)	120°C/min. (216°F/min.)	—
1450°C (2642°F) — 1600°C (2912°F)	10°C/min. (18°F/min.)	—
1600°C (2912°F)	—	20 min
1600°C (2912°F) — Room Temp. (*)	-120°C/min. (-216°F/min.)	—

* The restorations may be removed at 800°C (1472°F) or less depending on the circumstances.

IV. Composition

ZrO₂ (80-95%), Y₂O₃ (3-15%)

Units in parentheses are mass %.

V. Type and Class (ISO6872:2015)

Type:II/ Class:3

VI. Physical Properties

Coefficient of Thermal Expansion (25-500°C (77-932°F)): 9.7×10⁻⁶/K

VII. Directions for Use

- (1) Take the disc from the packaging and confirm that the disc does not have a crack or other damage.
- (2) Place the disc into the milling machine; then begin the milling process following the milling systems technical instructions.
- (3) After milling, remove the restorations from the disc with a diamond bur, etc.
- (4) The cutting waste or dust, which is attached to the restorations, can be removed with a gentle air stream.
- (5) Put the restorations into the refractory sagger tray and place them in the sintering furnace.
- (6) Depending on the performance of the sintering furnace used, review the sintering schedule shown above (III. Sintering Program) before sintering the restorations.
- (7) After sintering, adjust the restorations with a diamond bur as needed.
- (8) Confirm that the restorations have no cracks.
- (9)-1 FCZ:

Glaze baking: Create a high shine surface by polishing, especially on the contact areas, then apply the glaze on all surfaces in the usual manner.
- (9)-2 Frameworks: Build-up the dental porcelain (CERABIEN ZR or CZR PRESS LF) on the frameworks following the manufacturer's technical instructions. Check the coefficient of thermal expansion of the porcelain in the manufacturer's technical instructions to confirm compatibility.

VIII. Remarks on Handling

Contraindications:

1. If the patient is hypersensitive to zirconia or any other components, this product must not be used.
2. Do not use this product to make 4+ unit bridges, cantilever bridges or bridges containing molars.

Warning:

- If the patient or the dental professional demonstrates a hypersensitivity reaction, such as rash, dermatitis etc., discontinue use of the product and seek medical attention immediately.

Caution:

1. This product should not be used when malocclusion, clenching or bruxism conditions are present.
2. When milling the disc or cutting, grinding and polishing the restorations, use an approved dust mask and vacuum with air filter to protect your lungs from inhaling dust.
3. When milling the disc or cutting, grinding and polishing the restorations, use safety glasses to prevent dust from getting into your eyes. If the dust gets into your eyes, immediately rinse with copious amounts of water and consult a physician.
4. Do not use for any purposes except for dental restoration. This product is for dental application only.
5. The use of this product is restricted to dental professionals.
6. Do not touch the items heated in the furnace with your bare hands.
7. Do not take the restorations out of the sintering furnace during high temperature, as the quenching causes the breaking. However, if an automatic opening type furnace is used under the sintering program 2 or 3, the restorations may be removed from the furnace at 800°C (1472°F) or less. When removed from the furnace, the restorations must be put on a tray made of ceramic fiber (for example: Noritake Porcelain Mat) to cool slowly.
8. There may be a large difference between the setting temperature in the program and the temperature in the real furnace when using the furnace for sintering program 2 or 3. Please use the furnace after checking with the manufacturer that the furnace and temperature listed in sintering program 2 or 3 is compatible.
9. **Sintering program 2 or 3 is only recommended for frameworks (up to 3 units), FCZ crowns, FCZ bridges (up to 3 units), inlays, onlays and veneers.**
10. This product is particularly translucent; therefore, one should carefully consider its appropriateness for use when masking the underlying surface, such as a metal abutment, or a very dark or discolored tooth preparation.
11. The margins should be prepared with a deep chamfer and rounded shoulders, with cutting edges and corners rounded to eliminate sharp preparation corners. The angle of the axial surface should be within the range of 5 to 15 degrees.
12. When preparing teeth, avoid the following: deep shoulders, J-margins, knife edges, serrated margins, non-tapered abutments, undercuts, guide grooves, the formation of retentive holes, and sharp corners.
13. Keep the following thickness of this product for fabricating prosthetics:

Location & indication	Wall thickness
Anterior crown or bridge	0.8 mm or more
Veneer	0.4 mm or more*
Posterior crown or bridge	1.0 mm or more
Inlay or onlay	1.0 mm or more

* 0.4 mm or more of this product is for full zirconia veneers. Keep thickness 0.8 mm or more, if it is used for combination with the porcelain.

14. Use the following cross-sectional areas for connectors when fabricating bridges.

Location & indication	Connector cross section
Anterior 2- or 3-unit bridges	12 mm ² or more
Premolar 2- or 3-unit bridges	16 mm ² or more

15. Choose a shade color that is brighter than the intended color for a thick restoration, as it may look duller depending on the thickness of the restorations.
16. When using a sintering furnace for the first time and changing a sintering condition, colors after sintering may vary. Sinter a small piece of Zirconia beforehand and confirm the color.

Caution in conjunction with usage:

1. Do not use the disc if there is a crack noticed after removing it from the package.
2. Dry milling is recommended. If wet milling is used, the translucency of the disc may be reduced.
3. If you find a crack in the restorations, do not use.
4. The above sintering recommendation is only a guideline; some adjustments may be required depending on each individual furnace.
5. When milling the disc, use caution when approaching the milling bur with the plastic ring to prevent detachment of the disc. Mill the disc as if leaving a zirconia part contacting the internal side of the plastic ring at 2 mm.
6. Cut and remove the plastic ring prior to sintering the whole disc.
7. Dispose of this product as a medical waste to prevent infection.

Storage:

1. Store in a cool and dry place. Keep away from direct sunlight.
2. The product should be stored at 10-30°C (50-86°F).
3. Do not remove the disc from its packaging during storage.
4. The disc is fragile, and requires care when handling.
5. The product must be stored in an appropriate place where only dental personnel have access.
6. The product must be used by the expiration date indicated on the package.

[WARRANTY]


Kuraray Noritake Dental Inc. will replace any product that is proven to be defective. Kuraray Noritake Dental Inc. does not accept liability for any loss or damage, direct, consequential or special, arising out of the application or use of or the inability to use these products. Before using, the user shall determine the suitability of the products for the intended use and the user assumes all risk and liability whatsoever in connection therewith.

[NOTE]

If a serious incident attributable to this product occurs, report it to the manufacturer's authorized representative shown below and the regulatory authorities of the country in which the user/patient resides. "KATANA" and "CERABIEN" are registered trademarks or trademarks of NORITAKE CO., LIMITED.

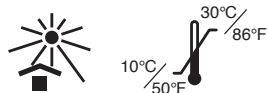UTML (Ultra Translucent Multi Layered)

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