

kuraray

Noritake

CLEARFIL™ Universal Bond Quick 2

The everything bond



Next-Generation
Single-Step Adhesive



Next-Generation Single-Step Adhesive

CLEARFIL™ Universal Bond Quick 2

CLEARFIL™ Universal Bond Quick 2 is the latest evolution in dental bonding technology, offering a streamlined, efficient solution for dental practitioners. Designed to enhance your workflow, this adhesive provides robust performance across a variety of clinical situations.



Choose CLEARFIL™ Universal Bond Quick 2

EASY APPLICATION

CLEARFIL™ Universal Bond Quick 2 is very easy to apply. The new filler dispersion method used makes it possible to enhance the products' flowability without lowering the filler content.

QUICK AND STRAIGHTFORWARD PROCEDURE

The streamlined procedure enables quick application without any waiting time or the need for rubbing. It involves minimal steps and is fast, making it well-suited for challenging cases where moisture control is difficult, or in pediatric dentistry.

CONSISTENTLY DURABLE BOND

Different well designed, matched ingredients create a solid bond to enamel and dentin.

HIGH STRENGTH BONDING LAYER

The added urethane tetra methacrylate monomer improves the strength of the bonding layer itself. The result is an increased resistance to shrinkage forces.

VERSATILITY FOR VARIOUS APPLICATIONS

Suitable for both direct and indirect restorations, as well as core build-up procedures, this adhesive provides the flexibility needed for a wide range of clinical scenarios. It is also compatible with self adhesive cements from other manufacturers, offering you greater choice.

MINIMIZED RISK OF POOLING

Enhanced flowability provides for even spreading, creating a thin, uniform bonding layer. This minimizes the risk of pooling thus maintaining the original morphology of the preparation.

READY TO USE ANYTIME

The adhesive can be stored at room temperature (2-25°C) and eliminates the need for shaking. Hence, it is always ready for immediate use, fitting seamlessly into your practice.




SCAN FOR MORE INFORMATION ABOUT CLEARFIL™ Universal Bond Quick 2



The heart of the universal portfolio

At Kuraray Noritake Dental Inc., we understand the demands of modern dental practices. Balancing complex procedures, tight schedules, and patient care is no easy task. That is why we have spent years refining our materials—so you can focus on what matters most: your patients. Our universal products simplify your workflow, offering you high-quality materials and tools that work seamlessly with your expertise.

CLEARFIL™ Universal Bond Quick 2 forms literally the connection between several universal resin-based products of our portfolio. Plus: there is no need for a separate dual-cure activator when used with other manufacturers' self-adhesive cements.

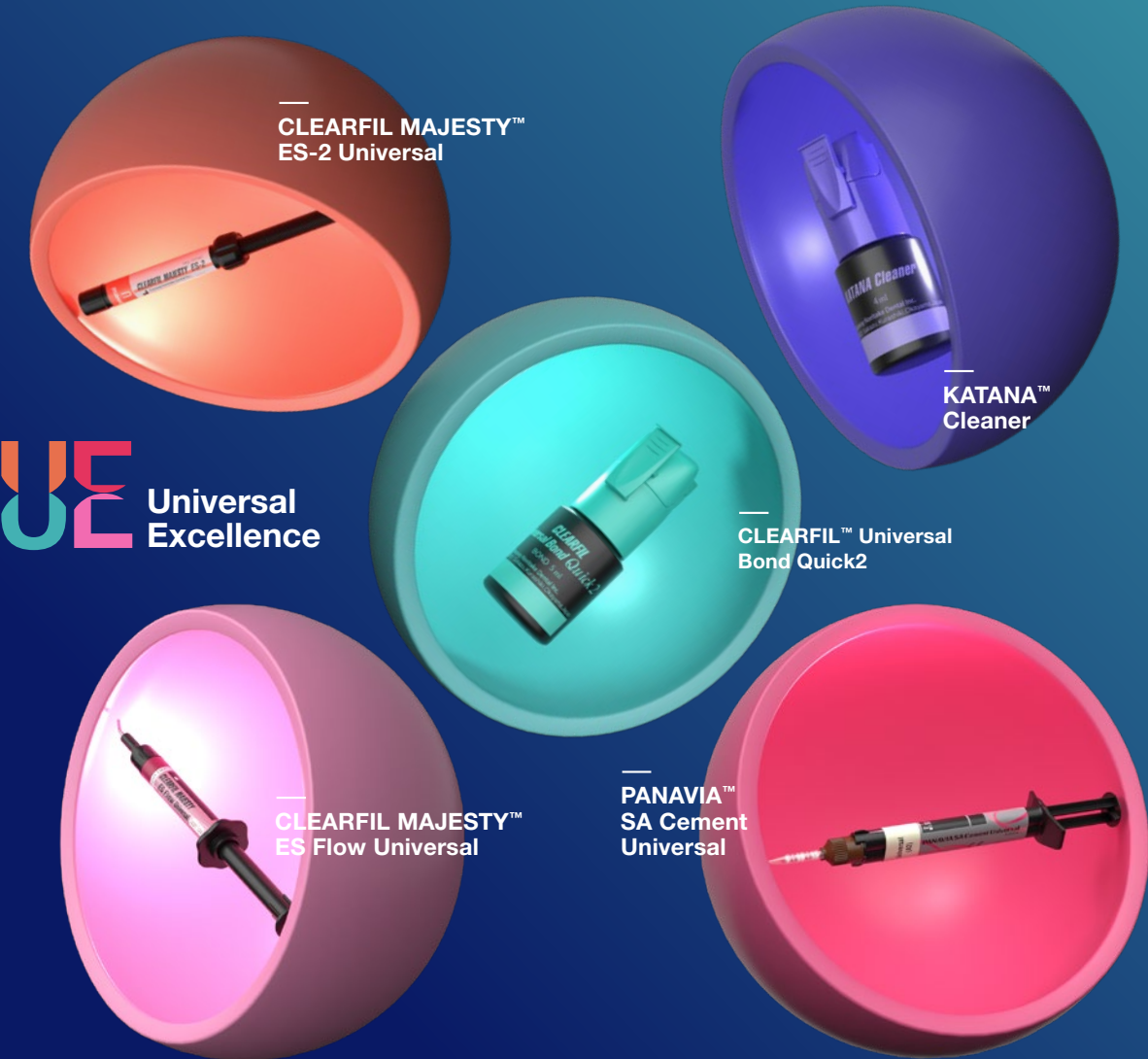

SCAN FOR MORE INFORMATION ABOUT Universal Excellence






The complex technology behind simplicity

At Kuraray Noritake Dental Inc., we have refined our monomer technology to deliver exceptional bonding performance with minimal effort. Our proprietary **Advanced Rapid Bond Technology** combines three key monomers to create a powerful adhesive that supports both ease of use and long-lasting results.



Together, these monomers form the foundation of **CLEARFIL™ Universal Bond Quick 2**, giving you a fast, dependable, and versatile adhesive for all your restorative needs.

Performance you can trust

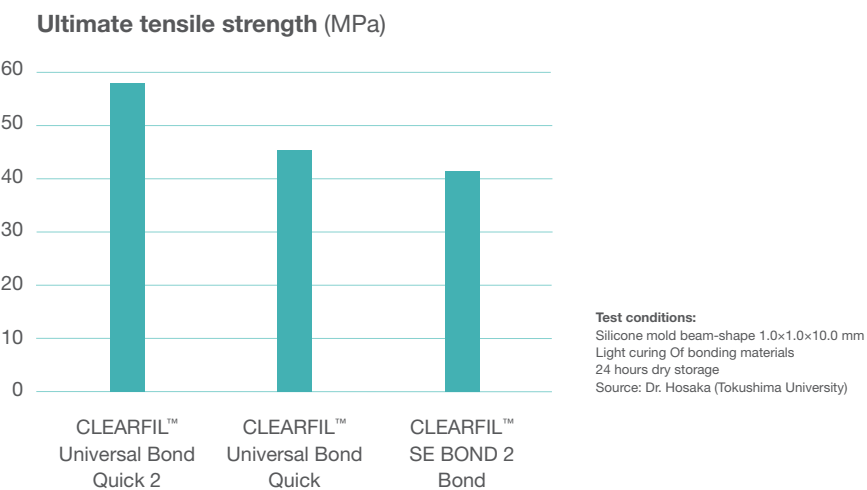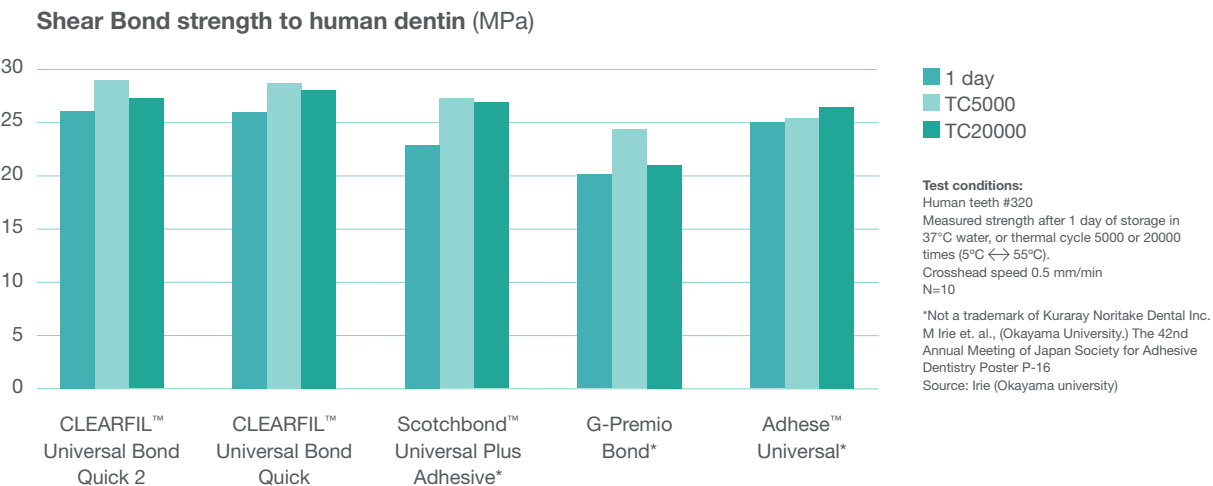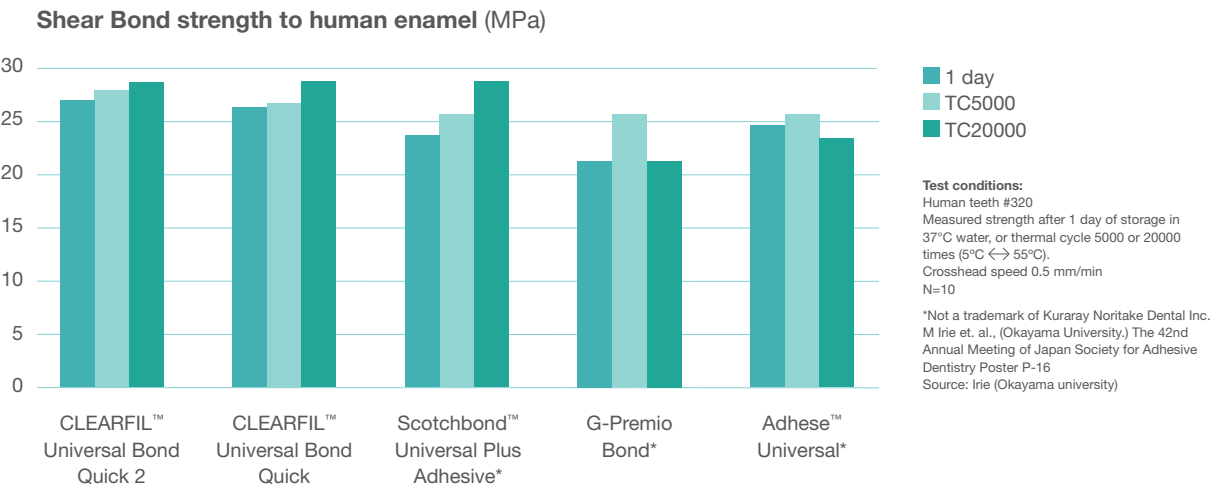

Direct restorations workflow

Follow the standard procedures for isolation, moisture control, cavity preparation and pulp protection

1 Tooth pretreatment Choose either etching procedure

a. Self-etching (Move to section 2)

b. Selective-etching Apply K-ETCHANT Syringe to the uncut and/or cut enamel, then rinse and dry

c. Total-etching Apply K-ETCHANT Syringe to the entire cavity (enamel and dentin), then rinse and dry

2 Apply BOND with a rubbing motion

3 Dry by blowing mild air until BOND does not move¹

4 Light-cure²

5 Place composite resin, light-cure and finish


Dental curing unit and curing time

Type	Light source	Light Intensity	Light-curing time
Halogen	Halogen lamp	More than 400 mW/cm ²	10 seconds
LED	Blue LED*	800-1400 mW/cm ²	10 seconds
		More than 1500 mW/cm ²	5 seconds

² Refer to table on the right for the light-curing time.

The effective wavelength range of each dental curing unit must be 400-515nm.
*Peak of emission spectrum: 450-480nm

How do we compare

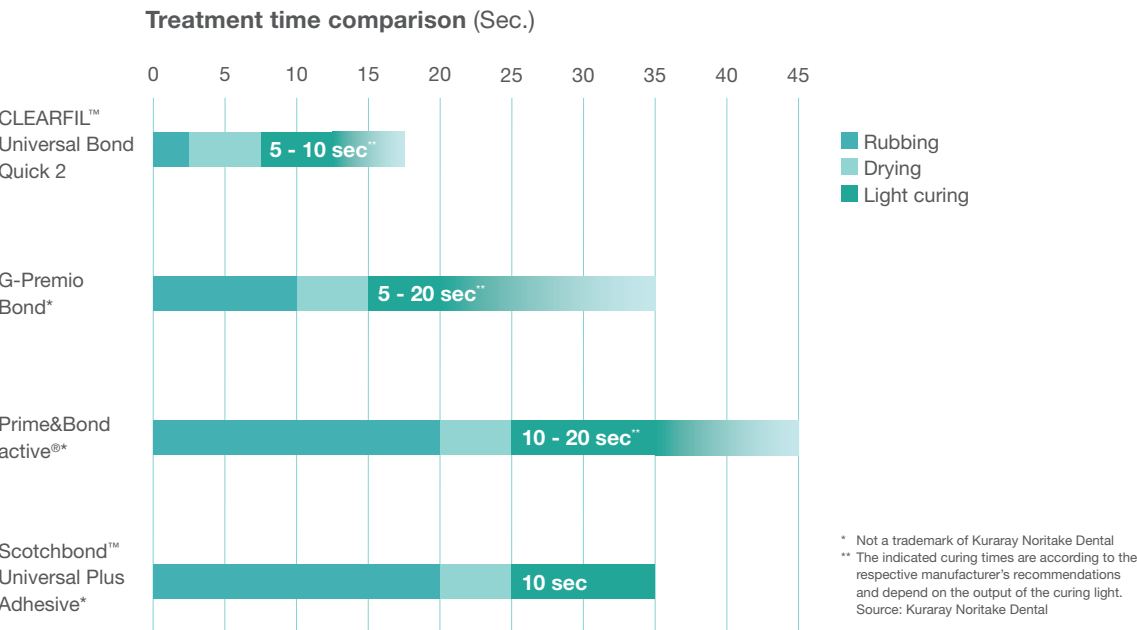

Ordering information

CLEARFIL™ Universal Bond Quick 2

#4042-EU	Bottle Refill	5.0 ml x 1 bottle
#4044-EU	Bottle Value Pack	5.0 ml x 3 bottles
#4046-EU	Unit Dose Standard Pack	0.1 ml tip x 50 tips, applic. brush fine silver x 50 pcs

UNIVERSAL EXCELLENCE PRODUCTS

CLEARFIL MAJESTY™ ES Flow Universal

#4600-EU	Low (U)	1 syringe 1.5 ml x 1, needle tip x 15 pcs
#4601-EU	Low (UD)	1 syringe 1.5 ml x 1, needle tip x 15 pcs
#4620-EU	Super Low (U)	1 syringe 1.5 ml x 1, needle tip x 15 pcs
#4621-EU	Super Low (UD)	1 syringe 1.5 ml x 1, needle tip x 15 pcs

CLEARFIL MAJESTY™ ES-2 Universal

#3181-EU	Universal (U)	1 syringe 2.0 ml
#3183-EU	Universal Light (UL)	1 syringe 2.0 ml
#3184-EU	Universal Dark (UD)	1 syringe 2.0 ml
#3182-EU	Universal White (UW)	1 syringe 2.0 ml
#3191-EU	Universal (U)	20 PLT 0.25g
#3193-EU	Universal Light (UL)	20 PLT 0.25g
#3194-EU	Universal Dark (UD)	20 PLT 0.25g
#3192-EU	Universal White (UW)	20 PLT 0.25g

PANAVIA™ SA Cement Universal

#4202-EU	Handmix Universal (A2)	1 syringe 5.0 ml
#4203-EU	Handmix White	1 syringe 5.0 ml
#4204-EU	Handmix Translucent	1 syringe 5.0 ml
#4212-EU	Automix Universal (A2)	1 syringe 8.2g, 20 Mixing tips
#4213-EU	Automix White	1 syringe 8.2g, 20 Mixing tips
#4214-EU	Automix Translucent	1 syringe 8.2g, 20 Mixing tips

KATANA™ Cleaner

#3970-EU	4.0 ml x 1 bottle
----------	-------------------

OTHER

K-ETCHANT Syringe

#3252-EU	2 Syringes 3.0 ml, needle tip x 40 pcs
----------	--



EU Importer

Kuraray Europe GmbH

Philipp-Reis-Strasse 4,
65795 Hattersheim am Main, Germany
Phone +49 (0)69 305 35 835
Fax +49 (0)69 305 98 35 835
www.kuraraynoritake.eu
centralmarketing@kuraray.com

- Before using this product, be sure to read the Instructions for Use supplied with the product.
- The specifications and appearance of the product are subject to change without notice.
- Printed color can be slightly different from actual color.

"CLEARFIL", "CLEARFIL MAJESTY" and "PANAVIA" are registered trademarks or trademarks of KURARAY CO., LTD.

"KATANA" is a registered trademark or trademark of NORITAKE CO., LIMITED.



Kuraray Noritake Dental Inc.
1621 Sakazu, Kurashiki, Okayama 710-0801, Japan
Website www.kuraraynoritake.com

