

# REALITY

RATINGS & REVIEWS | practical • unbiased • trusted  
36<sup>TH</sup> anniversary

## The Ratings | Desensitizers — Office

# Teethmate Desensitizer

### Manufacturer

Kuraray

<https://www.kuraraydental.com>

### Price

#### Set

\$135.06 (Since this is a powder/liquid, we cannot calculate its unit price/ml, but according to the manufacturer, each use costs \$0.90)

### Shelf life

2 years



### Raves & Rants

- ✓ Amazingly effective
- ✓ No harsh chemicals
- ✗ Powder-liquid formula requires hand-mixing
- ✗ Rinsing adds an extra step

## 2021

### REALITY

RATINGS & REVIEWS

Four Star Award



★★★★☆ 4.4

## Introduction/ Manufacturer's Claims

As long as our patients have sensitive teeth, manufacturers will continue to bring new desensitizers to the market. Kuraray, the company who has brought us such stalwart products such as Panavia and Clearfil SE Bond, has decided it was time to throw its hat into ring with Teethmate Desensitizer.

Unlike numerous other products in this category that have cloned the Gluma formulation of glutaraldehyde and HEMA, Teethmate Desensitizer takes a radically different approach. In the first place, it comes in two components: powder and liquid. Yes, that's correct, powder and liquid. And of course, you have to manually mix these components in a dispensing well.

As the story goes, when you mix the two components and apply it to teeth, it crystallizes to form hydroxyapatite, which presumably seals dentinal tubules and heals enamel fractures. It is being promoted virtually across the board, including for sensitive root surfaces, before/after bleaching, and under restorations. You don't have to evaporate any solvent, it doesn't require light curing, and it has no film thickness.

## Composition

**Powder** Tetra calcium phosphate (TTCP) and dicalcium phosphate anhydrate (DCPA).

**Liquid** Water and a preservative.

Most (65%) evaluators liked the fact that it does not contain strong chemicals, while the other 35% were more concerned with its effectiveness regardless of what it contains. One evaluator noted that it is a truly biomimetic approach to treating sensitivity and questioned why would you want to use polymers or other agents when you can use the "real" thing, aka hydroxyapatite? On the other hand, another evaluator stated that using chemicals like glutaraldehyde and HEMA haven't caused much of a problem, so why the concern?

## pH

10 (according to the manufacturer).

## Viscosity

The thick viscosity mixture was deemed to be acceptable by most (90%) evaluators, with only 10% considering it to be too thick.

## Odor

Most (65%) evaluators thought it has virtually no odor, while the other 35% found it to be about the same as other desensitizers.

## Color

White.

## Use

Clean the tooth as you would with any product. All evaluators except one blotted the sensitive areas after rinsing the cleaner, leaving them moist, while the lone holdout dried the teeth.

Then, for one tooth, you dispense one small scoop of powder into the supplied dispensing well and add one drop of liquid. For 2-3 teeth, use the large scoop and two drops of liquid. With the supplied applicator tip (very small), you mix the components for at least 15 seconds to produce a relatively thick paste, which you then apply to the teeth for 30 seconds using a rubbing-like motion, followed by rinsing it off. If necessary, you can repeat this sequence several times.

Mixing is surprisingly easy – the powder wets very quickly and producing the slurry is almost automatic. Applying it to the tooth is also easy. Since it's a paste, it's not runny and rubbing it into the tooth is not an issue. And even though the paste is white, rinsing it off basically leaves the tooth without any residual material, thus confirming the claim of no film thickness.

As far as dispensing powder and liquid, most (60%) evaluators did not consider it a big inconvenience, while the other 40% thought it was cumbersome. Mixing wasn't a major issue for 75% of the evaluators, while the other 25% considered it to be cumbersome. One evaluator noted that mixing with the applicator tip was not very efficient and you needed to be sure the bottle of liquid was perfectly vertical to get the drops in the correct volume.

Rubbing/scrubbing it into the sensitive area was the preferred application method of most (65%) evaluators, while the other 35% used gentle agitation. Although 30 seconds of application was most typically used (75%), a few evaluators applied it longer and several used more than one application.

And, despite the thickness, most (80%) evaluators thought applying it was easy, especially since it isn't runny. Nevertheless, the other 20% still preferred a liquid desensitizer.

Rinsing it off was also found to be easy by most (90%) evaluators, while the other 10% just thought the rinsing step required too much effort.

## Effect on Bond Strength

Slight positive effect.

## Effectiveness

### **Under restorations placed with etch/rinse adhesives**

Most (75%) evaluators reported it worked all the time, while the other 25% found it worked some of the time. One evaluator noted that this parameter is more difficult to verify since he rarely gets sensitivity under restorations anyway.

## Cervical Areas

Half of the evaluators reported it worked all the time after 2-3 applications, while another 15% found it worked all the time after only one application. The other 35% found it worked some of the time. One evaluator noted

that patients noticed a difference only after the second application.

## Before Bleaching

This was a mixed bag, with half of the evaluators reported it worked all the time, 17% some of the time, and 33% did not find it to be effective at all.

## After Bleaching

Most (67%) evaluators reported it worked all the time, while the other 33% found it worked some of the time.

## Packaging

The plastic shrink-wrapped kit comes in a conventional cardboard box with product identification on the top and three sides, with the expiration date on one side. The powder is in a typical small glass jar with a scoop leveler, while the liquid is in a conventional plastic squeeze bottle.

Virtually all evaluators thought the packaging was nice, neat, and organized with good access to the components. However, two evaluators felt keeping the contents aseptic left something to be desired.

## Directions

Annoying plain paper foldout design in three languages and nano-type. There is also a plastic-coated treatment room card. The paper instructions are somewhat redundant but fairly easy to follow. The treatment room card features color illustrations and easy-to-follow, step-by-step application procedures. It is well done. However, one evaluator felt the cleaning instructions could have been more detailed.





## Strengths

Once it is mixed, it is easy to apply and rinse. Effective most of the time. Mixture is easy to see and control since it is white and not runny. Chemistry is very biocompatible and biomimetic — no harsh chemicals, which can be especially important in patients with chemical sensitivities. Has slight positive effect on bond strength with etch & rinse adhesive. Minimal if any odor. No film thickness..



## Weaknesses

Doesn't always work. Dispensing and hand-mixing a powder-liquid material is a nuisance. Rinsing it off is an additional step. Packaging has asepsis challenges.

## BOTTOM LINE

**Although hand-mixing a powder-liquid product is definitely not as easy as dispensing a premixed, single-component liquid or gel, its common sense approach of recreating hydroxyapatite with the sophisticated chemistry that does not employ a controversial ingredient such as glutaraldehyde definitely has appeal and its effectiveness is impressive.**



**Michael B. Miller, D.D.S. – President/Editor-in-Chief**

**Ingrid R. Castellanos, C.D. – Vice President/Publisher**

*Editorial Team:* M. Bilge Akbulut, D.D.S., Ph.D., Konya, Turkey; Leslie C. H. Ang, B.D.S., M.Sc.; Robert W. Baker, Jr., D.M.D, Ithaca, NY; Sema Belli, D.D.S., Ph.D., Konya, Turkey; Nathan S. Birnbaum, D.D.S., Wellesley, MA; Alan A. Boghosian, D.D.S., Chicago, IL; Sonia Regina Bordin-Aykroyd, D.D.S., M.Sc., Sao Paulo, Brazil; Matthew Brock, D.D.S., C.A.G.S., M.S.D., Chattanooga, TN; Richard L. Brown, D.D.S., Tulsa, OK; Mitch A. Conditt, D.D.S., Fort Worth, TX; Simona Cuevas, D.D.S., San Antonio, TX; Marvin A. Fier, D.D.S., Pomona, NY; George A. Freedman, D.D.S., Toronto, Ont., Canada; John Gammichia, D.M.D., Orlando, FL; Fay Goldstep, D.D.S.; Toronto, Ont., Canada; Gary Henkel, D.D.S., Horsham, PA; David S. Hornbrook, D.D.S., San Diego, CA; Timothy F. Kosinski, M.S., D.D.S., Bingham Farms, MI; Hannu O. Laamanen, D.D.S., M.S., Turku, Finland; Carole Landman, D.D.S., Chicago, IL; Clarence C. Lindquist, D.D.S., Washington, D.C.; Edward Lynch, M.A., B.D.Sc., Ph.D., London, UK; Hans Malmstrom, D.D.S., Rochester, NY; Robert C. Margeas, D.D.S., Des Moines, IA; Sandesh Mayekar, M.D.S., Mumbai, India; Steven McGowan, C.D.T., Seattle, WA; Michael K. McGuire, D.D.S., Houston, TX; Elaine Mo, B.D.S., London, UK; Aikaterini Papathanasiou, D.D.S., Boston, MA; Christopher Pescatore, D.M.D., Danville, CA; Stephen D. Poss, D.D.S., Brentwood, TN; Robert G. Ritter, D.M.D., Juniper, FL; Joseph Sabbagh, D.D.S., M.Sc., Ph.D, Beirut, Lebanon; Andrew T. Shannon, D.D.S., Vancouver, BC, Canada; Todd Snyder, D.D.S., Laguna Nigel, CA; Liviu Steier, D.M.D., Needham, MA; Franklin Tay, B.D.Sc.(Hons), Ph.D., Augusta, GA; Stephanie Tilley, D.M.D., Pensacola, FL; Brian Toogood, D.D.S., Abilene, TX; Marcos A. Vargas, D.D.S., M.S., Iowa City, IA; Nicholas Vece, D.M.D., Mt. Kisco, NY; Charles Wakefield, D.D.S., Lewisville, TX; Thomas G. Wilson, Jr., D.D.S., Dallas, TX; David Winkler, D.D.S., Odense, Denmark; Tyler Wynne, D.D.S., Clemmons, NC.

**To become a member of REALITY RATINGS & REVIEWS, please visit our Web site at [www.realityratings.com](http://www.realityratings.com).**

#### **NO COMMERCIALIZATION POLICY**

We accept no advertising and are not beholden to any commercial interest. Product evaluations and ratings are intended only to guide our readers to make wise and informed purchases. The unauthorized use of product evaluations and ratings in advertising or for any other commercial purpose is strictly forbidden.

**REALITY RATINGS & REVIEWS** (ISSN#1041-8253) is an online information service from **REALITY** Publishing Company, 1322 Eagle Point Dr., Georgetown, TX 78628, U.S.A., 512-688-3423, Fax 512-489-9225. A one-year membership includes access to the online database plus periodic PDF issues of **REALITY NOW**. Call for membership and publication rates or access our Web site for enrollment information. Payments by check must be in U.S. funds drawn on a U.S. bank, or by Visa, MasterCard, or American Express. All rights reserved. No part of **REALITY RATINGS & REVIEWS** or **REALITY NOW** may be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopying, recording, or by any information storage and retrieval system without the written permission of the Publisher, except where permitted by law. Copyright ©2021 by **REALITY** Publishing Company. GST #898-896-659. POSTMASTER: Send address changes to **REALITY** Publishing Company, 1322 Eagle Point Dr., Georgetown, TX 78628.



A publication member of the American Association of Dental Editors

**kuraray**

## TEETHMATE™ DESENSITIZER

Creating  
Hydroxyapatite

That's How You Treat Sensitivity



**KURARAY AMERICA, INC.**

33 Maiden Ln, Suite 600D, New York, NY 10038

800 879-1676 | [kuraraydental.com](http://kuraraydental.com)

©2013 Kuraray America, Inc.



## OPTIMUM TREATMENT OPTION

**TEETHMATE™ DESENSITIZER** (TMD) consists of a powder and liquid. We optimized delivery by using powder and liquid which are hydrophilic components. The product readily penetrates into the moist tooth structure

to occlude opened dentinal tubules or enamel cracks by forming hydroxyapatite (HAp) and usually does not require multiple treatments to achieve a therapeutic effect.




## EASY TO MIX AND APPLY TO ANY SURFACE


Application is easy. Clean the surface, saliva or plaque; hold the liquid container vertically and dispense the liquid (Figure 6);

add the powder in the correct ratio (Table 1) and mix with a micro-applicator. Rub the slurry into the treatment area for at

least 30s. Rinse and check for desensitizing effect.

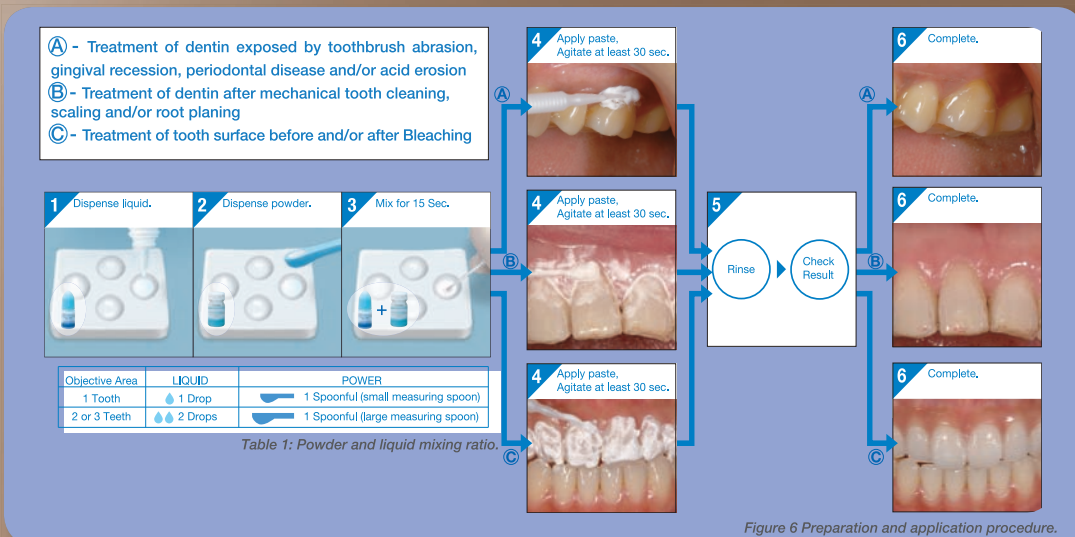

## PRODUCT INFORMATION

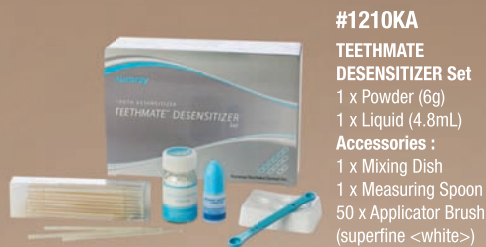

**#1210KA**  
**TEETHMATE**  
**DESENSITIZER Set**  
 1 x Powder (6g)  
 1 x Liquid (4.8mL)  
**Accessories :**  
 1 x Mixing Dish  
 1 x Measuring Spoon  
 50 x Applicator Brush  
 (superfine <white>)



**#1216KA**  
**TEETHMATE**  
**DESENSITIZER**  
**POWDER**  
 (6g)



**#1217KA**  
**TEETHMATE**  
**DESENSITIZER**  
**LIQUID**  
 (4.8mL)

Warranty: Kuraray Noritake Dental Inc. will replace any product that is proved to be defective. Neither Kuraray Noritake Dental Inc. nor the distributor accepts liability for any loss or damage, direct, consequential or special, arising out of the application or use of or the inability to use these products. Before using, the user should determine the suitability of the products for the intended use and the user assumes all risk and liability whatsoever in connection therein. TEETHMATE DESENSITIZER and Kuraray Noritake Dental Inc. products are available only through Kuraray Authorized dealers.

"CLEARFIL" and "TEETHMATE" are trademarks of KURARAY Co.,Ltd.

Manufactured by Kuraray Noritake Dental Inc.  
 1621 Sakazu, Kurashiki, Okayama 710-0801, Japan

1210LT