

**kuraray**

*Noritake*

**ALL THAT REMAINS  
IS A GOOD FEELING**

**CLEARFIL™ SE Protect**






## YOUR NEW ADHESIVE WITH AN ANTIBACTERIAL EFFECT

CLEARFIL™ SE Protect is the next version of the world's first self-etching adhesive with a reliable antibacterial effect. CLEARFIL™ SE Protect contains the functional monomer MDPB, which exhibits an "Antibacterial Cavity Cleansing Effect."

With its MDPB monomer, CLEARFIL™ SE Protect destroys bacteria on contact. How does this bacterial destruction work? Quite simply: the cell membrane of bacteria is usually negatively charged, the MDPB monomer is positively charged, hence drawn to the negative contact point of the cell membrane. At this point, the bacteria lose their electrical balance and the cell membrane of the bacteria is destroyed, by literally "bursting". This process, also known as bacteriolysis, kills the bacteria.

CLEARFIL™ SE Protect also contains Kuraray Noritake Dental's unique monomer MDP, which strongly bonds to hydroxyapatite. Thanks to these MDPB and MDP monomers, CLEARFIL™ SE Protect leads to long-lasting and reliable restorations.

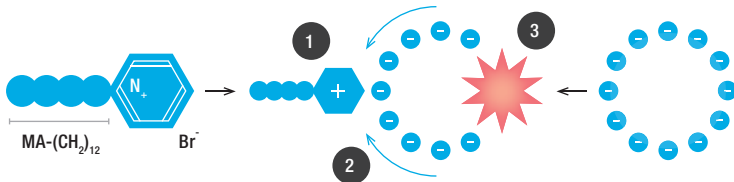

CLEARFIL™ SE Protect



CLEARFIL™ DC Activator

## ORDER INFO

### CLEARFIL™ SE Protect

**Kit #2870-EU:** Primer (6 ml), Bond (5 ml), 50 x 2 disposable brush tips (white), 1 brush tip handle (white), 1 brush tip handle (black), 1 mixing dish, 1 light blocking plate, 1 outer case;

**Primer #2882-EU:** 1 bottle (6 ml);

**Bond #2881-EU:** 1 bottle (5 ml);

**Value Pack #2872-EU:** 3 Primer (6 ml), 3 Bond (5 ml);

## PRODUCT RECOMMENDATION

### CLEARFIL™ DC Activator

**#3250-EU:** 1 bottle (4ml);

## YOUR CONTACT

Kuraray Europe GmbH  
BU Medical Products  
Philipp-Reis-Strasse 4  
65795 Hattersheim am Main

**Phone** +49 (0) 69-305 35 833  
**Fax** +49 (0) 69-305 98 35 833  
**E-Mail** dental@kuraray.eu  
**Website** www.kuraraynoritake.eu